



Portable Oxygen Concentrator

Standing at just 49cm tall, the SeQual® Eclipse 5™ is one of the most technologically advanced oxygen concentrator systems in the world today.

The on-board control panel is feature rich, yet easy to use offering both continuous flow and pulse dosing options, great for use during the day and whilst sleeping at night.



Pulse Flow

1 to 9
Pulse flow
settings



Continuous

0.5 to 3
Litres Per
Minute



Weight

18.4lbs
(8.3kgs)



Battery life

Up to
5
Hours
on setting 2
pulse flow*

intermedical
H O M E C A R E

www.portableoxygen.co.uk | 0800 028 2194

Trusted in Respiratory Health

The all-in-one portable oxygen solution with advanced features.

The Eclipse 5™ comes with a “fine-tuning option” tailoring the speed of which oxygen is delivered on pulse dose in accordance to the patient’s respiratory requirements. Coupled with this is the Eclipse’s autoSAT[®] technology; a pulse delivery system which maintains a consistent pulse dose volume for every breath as the patient’s respiratory rate changes.

With the option to toggle between pulse flow and continuous flow, patients who suffer from sleep apnoea will certainly benefit from this all-in-one solution.

With its robust design and an array of optional accessories available, the Eclipse 5™ is the ideal travel companion for users preferring an all-round solution.




◀ Easy to use control panel with LCD display showing battery remaining and the current setting. Settings can be changed easily between pulse and continuous modes as well changing the flow by using the plus and minus buttons.



◀ The battery for Eclipse 5™ slots neatly into the base of the machine. The battery can last up to 5.4 hours on setting 1 (pulse flow mode) and it can be re-charged whilst in the machine or independently using the desktop charger (as shown - optional accessory).

Standard Package


- SeQual Eclipse unit
- Power cartridge
- Cart with telescopic handle
- AC mains power supply unit
- DC power supply unit
- Accessory bag
- User manual
- Nasal cannula

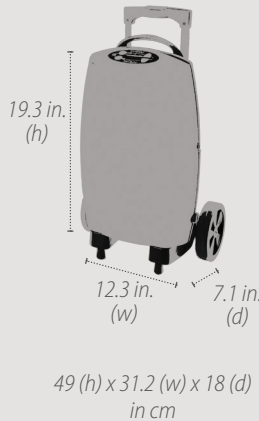

Approximate battery run-time
(in hours:minutes)

Continuous Flow Mode		Pulse Flow Mode	
0.5 lpm	4:25	Setting 1	5:25
1 lpm	3:40	Setting 2	5:05
2 lpm	2:00	Setting 3	4:55
3 lpm	1:20	Setting 4	4:00
		Setting 5	3:40
		Setting 6	3:30
		Setting 7	2:30
		Setting 8	2:00
		Setting 9	1:40

*Run-time on pulse flow mode dependent on breathing rate

Concentration	90% -3% /+ 5.5%
Weight	8.3 kg (18.4lbs) with power cartridge installed
Power	AC: 100-240v DC: 12v
Noise	40db(A) on setting 3 (pulse)
Batt. charge time	Approx. 2 - 5 hours
Warranty	3 years 1 year on battery & accessories
 Airline approved	Yes, subject to airline's terms & conditions



This device may be suitable for use whilst sleeping. This device may be used with a CPAP or BiPAP device.

intermedical
HOMECARE

0800 028 2194
www.portableoxygen.co.uk

Intermedical (UK) Limited, Homecare Division
Unit 6 Mill Hall Business Estate, Aylesford, Kent, ME20 7JZ, United Kingdom



Cert: 87M18032